

AMT MiniONE®

Capsule Non-Balloon & Non-Balloon Devices



Innovative Technology

The AMT MiniONE® Capsule Non-Balloon Button combines the advantages of a low-profile device with the longer life of a non-balloon device.

- Replacement occurs less often with no balloon to fail
- Capsule material is self-lubricating which facilitates ease and comfort on insertion
- The pull-away tether releases outer capsule to easily deploy internal bolster once placed

Improve Comfort

Comfort is key when it comes to patient, caregiver and clinician satisfaction with an enteral feeding device.

- Capsule design can reduce discomfort during insertion
- Comparable button shaft to a balloon button
- Ratcheted removal tool elongates bolster for removal
- Internal bolster is smaller than a balloon device
(less likely to block the pylorus)




- Secure Cap
- Glow Green™ Feed Port
- Soft External Bolster
- Tapered Transition
- Bolster (Post Deployment)



Pull for Capsule Deployment

Easier Nighttime Connections

The MiniONE® Capsule Non-Balloon and MiniONE® Non-Balloon feed port glow in the dark with AMT's **Glow Green™** Technology!

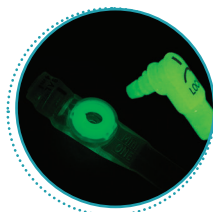

Innovating. Educating. Changing Lives.™

AMT MiniONE®

Capsule Non-Balloon & Non-Balloon Devices







Capsule Non-Balloon Size Availability

Diameter	Stoma Length (cm)								
14F	1.0	1.2	1.5	1.7	2.0	2.5	3.0		
18F		1.2	1.5	1.7	2.0	2.4	3.0	3.4	4.4
20F				1.7		2.4		3.4	4.4
24F				1.7		2.4		3.4	4.4

Non-Balloon Size Availability

Diameter	Stoma Length (cm)								
14F	1.0	1.2	1.5	1.7	2.0	2.5	3.0		
18F		1.2	1.5	1.7	2.0	2.4	3.0	3.4	4.4
20F				1.7		2.4		3.4	4.4
24F				1.7		2.4		3.4	4.4

Compatible Feed Sets

		Connector Options					
		Right Angle (GG)			Straight (GG)		
		2"	12"	24"	2"	12"	24"
Adapter Options	 Purple Y-Port		8-1255 8-1257G*	8-2455			
	 Bolus		8-1221	8-2421	8-0211	8-1211	8-2411
	 Clear Y-Port	8-0222	8-1222	8-2422	8-0212	8-1212	
	 Hybrid Y-Port		8-1255-H	8-2455-H			
	 ENFit® Y-Port		8-1255-ISOSAF	8-2455-ISOSAF			
	 ENFit® Single Port		8-1221-ISOSAF	8-2421-ISOSAF	8-0211-ISOSAF	8-1211-ISOSAF	8-2411-ISOSAF

GG=Glow Green™

■ Included in MiniONE® Legacy Kits

■ Included in MiniONE® ENFit® Kits

* Y-Port with Enhanced Med Port Plug

All come in a package of 10/box.

All MiniONE® feed sets are compatible with Mic-Key® devices. AMT will continue to offer our legacy product line along with a complete ENFit® product line. Complete MiniONE® ENFit® kit part numbers will end with a "-1", as a suffix for ENFit®, while legacy part numbers will not change. **Note:** Feed sets are supplied assembled. Choose the best combination of length, connector and adapter type for your use. Contact your insurance provider regarding the number of feed/extension sets allowed per month. LEGACY MiniONE® KITS are packaged with one 8-1255 Y-Port feed set and one 8-1211 Bolus feed set. ENFit® MiniONE® KITS are packaged with one 8-1255 Y-Port feed set and one 8-1211 Bolus feed set, along with two ENFit® syringes.

Spare Removal Tool

(For Obturated Non-Balloon & Capsule Non-Balloon Buttons)

AMT Obturator for elongating button during insertion and removal



Order #	Description	Qty/Box
R-S1012	14/1.0, 14/1.2	1
R-S1517	14/1.5, 14/1.7	1
R-S2023	14/2.0	1
R-S2427	14/2.5	1
R-S3033	14/3.0	1
R-L1012	18/1.2	1
R-L1517	18/1.5, 18/1.7, 20/1.7, 24/1.7	1
R-L2023	18/2.0	1
R-L2427	18/2.4, 20/2.4, 24/2.4	1
R-L3033	18/3.0	1
R-L3437	18/3.4, 20/3.4, 24/3.4	1
R-L4044	18/4.4, 20/4.4, 24/4.4	1

800 869 7382

CS@AppliedMedical.net



Applied Medical Technology, Inc.
8006 Katherine Boulevard
Brecksville, OH 44141
800 869 7382

www.AppliedMedical.net

An ISO 13485 Registered Company

©2018 Applied Medical Technology, Inc. C-4550-D

Patents: <http://www.AppliedMedical.net/Patents/> - Trademarks: <http://www.AppliedMedical.net/Trademarks/>

Mic-key is a registered trademark of Avent, Inc.